KariesScreenTest („KST“) MD2426; KariesScreenTest +P („KST +P“) MD2427. Instructions for Use (IFU; English language)

Introduction

The KST *after Laurisch et al.* is used for the identification and semiquantitative count of *streptococcus mutans* and *lactobacillus species* from human saliva specimen applying selective culture media *in vitro*.

A modified Rogosa-Agar (light color) is used to select for *lactobacillus species*.

A modified Mitis Salivarius-Agar (blue color) is used to select for *streptococcus mutans*.

Indication

KST is used to identify bacteria species that are strongly associated with caries development. KST is an *in vitro* diagnostic to be used by skilled medical personnel, only.

Warnings

For medical professional use only.

Store out of the reach of children.

Read these instructions for use and the safety data sheets before using KST.

Contraindications

Patient sample for the performance of KST may not be taken (1) during local/oral or systemic antibiotic therapy and two weeks thereafter, and not (2) within 12 hours after an antibiotic mouth rinse.

Patients must not do one of the following things within one hour prior to the specimen collection: -no drinking, no eating; -no chewing gum; -no smoking; -no teeth brushing.

Test Procedure

Take out one vial per patient from the box in the fridge and check the agar surfaces excluding contamination or bacterial growth. Warm up to RT before use (30 min. at RT or 5 min at 37°C).

KST should be performed by skilled medical personnel, only, in accordance with these instructions for use and good medical practice.

(1) Sample collection: *For patients 6 years and older*, let the patient chew one of the enclosed paraffin pellets over 5 minutes to stimulate salivation. Collect the saliva in a suitable container. *For patients 5 years or younger*, collect saliva and/or plaque specimen using a fine brush or swab.

(2) Remove the agar carrier from the test vial.

(3) Place a NaHCO₃-tablet at the bottom of the vial to generate appropriate culture conditions. If test is performed with plaque only or very little saliva, add two drops tap water onto the tablet.

(4) Carefully remove both foils protecting the growth media without touching the agar.

(5) Thoroughly, inoculate both agar media with saliva / the collected specimen wetting the entire agar surfaces. Do not scratch or punch the media surface.

(6) Slide the agar carrier back into the vial and screw-close the vial tightly.

(7) Mark the vial or the lid of the agar carrier with the patient name / code ensuring unique patient identification.

(8) Incubate the vial upright (lid up) over 48h at 36 ± 1°C (approx. 99°F). For appropriate incubators, please refer to www.aurosan-shop.de or contact us for advise.

(9) Immediately after removing the vial from the incubator analyze the bacterial growth by comparing the density of the colony forming units (CFU) with the respective reference pictures below.

Interpretation of results

Presence of both streptococcus mutans and lactobacillus in high numbers (>10⁵ CFU) are associated with a high risk to develop caries.

Therapeutic and prophylactic measures should be discussed between dentist and patient. A communication aid that can be individualized is available at www.revois-dental.com

Test performance and quality control

The growth media of KST are highly selective. The test may be interpreted only against the respective reference picture. For internal quality assurance the following control strains may be used: *Streptococcus mutans* ATCC™ 25175 (e.g. AUROSAN order n° MD0266P); *Lactobacillus acidophilus* ATCC™ 4356 (e.g. AUROSAN order n° MD0243P). For the performance of quality controls, the necessary equipment and training, please refer to the AUROSAN mediathek at www.aurosan.de or contact us, directly.

Test disposal

Single use only! Contaminated or used vials need to be disposed as biohazards in accordance with applicable national or local laws. Some regulations may permit the disposal in regular trash after autoclavation. Appropriate autoclaving bags are available (e.g. AUROSAN order n° MD-AL-00501).

Storage conditions

KST should be stored at 2-8°C (36-46°F) – do not freeze. Protect from light and humidity. Do not use after the expiration date printed on the agar carrier.

Package Content and Order Information

Manufacturer: AUROSAN GmbH, Frankenstr. 231, D-45134 Essen, Germany
Tel +49 (0)201 21961 701
Fax +49 (0)201 21961 731
E-Mail KariesScreenTest@aurosan.de
www.aurosan.de | www.revois-dental.com

Content MD2426 – KariesScreenTest:
10 vials with 2 media, 10 NaHCO₃ tables, 10 paraffin pellets, IFU
NaHCO₃ tables and paraffin pellets can be purchased separately, as well.

Content MD2427 – KariesScreenTest +P:
same as MD2426; in addition 10 buffering test strips.

Extensive product information online available at www.kariesscreentest.de
www.revois-dental.com
www.revois-shop.de

Order information for Dentists in Germany and Austria – please refer to www.minilu.de
Order information for Pediatricians – RÖSCH Medizintechnik - www.r-mt.de

Liability disclaimer

This *in-vitro* diagnostic is to be used by skilled medical personnel, only, and within the indication designated in this IFU. Processing, application and interpretation should strictly follow this IFU.

Therapeutic or prophylactic measures recommended or initiated by the dentist as an interpretation of KST are the sole responsibility of the dentist. The user of KST is solely responsible for the testing the material for its suitability and use for any purpose not explicitly stated in this IFU. Description, explanation, data or citation do not constitute any warranty of attributes and are not binding.

KST is made in Germany.

AUROSAN is certified in accordance with ISO13485.

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

MD2426 and MD2427: Stimulating salivation and collecting saliva: Have patient chew on paraffine pellet over 1 min. and swallow the saliva. Have him then chew for another 5 min. and collect saliva in a cup.


With a pipette or syringe, determine the volume of saliva per min. ((normal ≥ 1 ml/min; consider hyposalivation at values $\leq 0,5$ ml/min).

MD2427: Using a pipette, transfer a drop of the saliva onto the buffer test strip. (pls. also refer to backpage)

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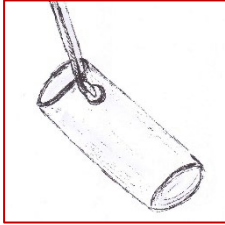
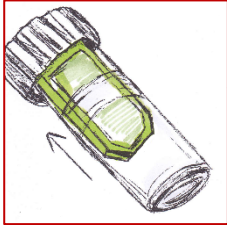
Remove the agar carrier from the tube and add a NaHCO₃-tablet into the tube. Remove protection foils from both sides without touching or damaging the agar surfaces.

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
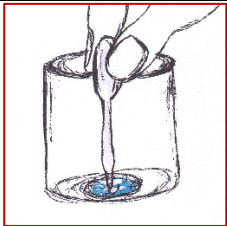
Using a pipette or syringe, transfer the collected saliva and entirely wetten both agar surfaces. Let excessive saliva drop off.

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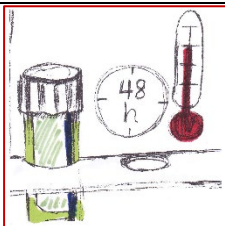
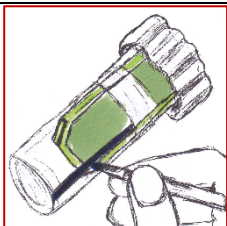
Put agar carrier back into the tube and thoroughly twist-close the tube. Don't forget to lable tube or top with patient ID. Inkubate for 48h at 36 ± 1 °C.

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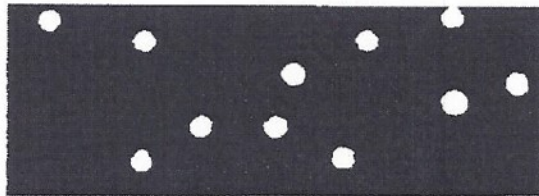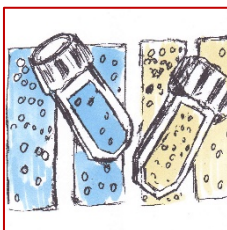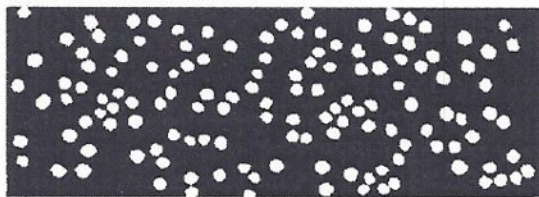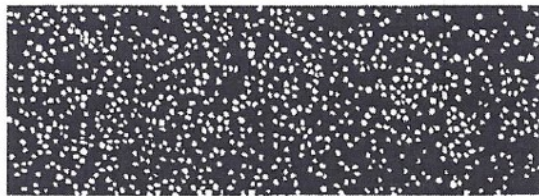
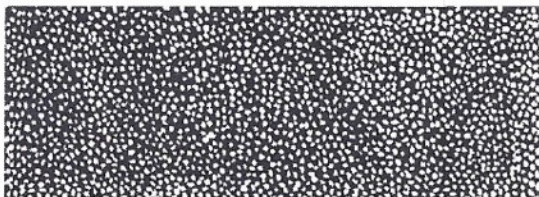
For the orientating measurement of colonies` density compare growth on agar media with reference pictures on the left: Document the result. Streptococcus colonies may be very small and shiny – consider using a magnifier.




Densities of 10⁵ CFU/ml (CFU = KBE) or more are correlated with a high risk to develop caries. A lowering of the saliva pH value or a reduced buffering capacity of the saliva (MD2427 only) points to an even higher risk.

KariesScreenTest - MD2426 / KariesScreenTest +P - MD2427
 Kurzanleitung für Patienten ab 6 Jahren – *short protocol for patients 6 years and older*

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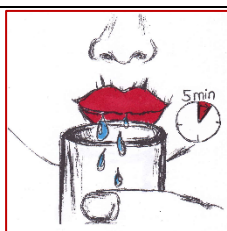*MD2426 and MD2427: Stimulating salivation and collecting saliva: Have patient chew on paraffine pellet over 1 min. and swallow the saliva. Have him then chew for another 5 min. and collect saliva in a cup.
 With a pipette or syringe, determine the volume of saliva per min. ((normal ≥ 1 ml/min; consider hyposalivation at values $\leq 0,5$ ml/min).*

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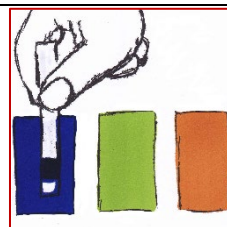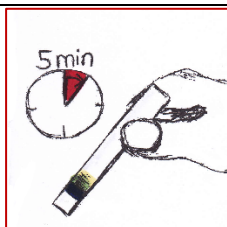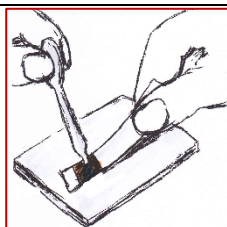*For MD2427 only, in addition:
 Using a pipette or similar, transfer a drop of the saliva onto the yellowish test field of the buffer strip.
 Wait 5 min.
 Compare color of the test field with the reference color pictures on the left. Blue color proves high, yellowish color low buffering capacity of the saliva.
 With proven low buffering capacity, explore additional caries preventive measures.*

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KariesScreenTest wird hergestellt von AUROSAN GmbH in Deutschland.
KariesScreenTest is manufactured by AUROSAN GmbH and produced in Germany.